
TitraLab® AT1000 and KF1000 Series Workstations

Software Update Procedure

1. On the PC, copy the content of 'USB Drive' folder from the zip file onto a USB flash drive.
The HACH\SOFT folder must contain a [Model name]_Vx_xx.CAB file.



2. Make sure that there are no USB flash drives connected to the instrument. Then, connect the USB flash drive containing the new software to the instrument.



3. Power on the instrument and press the down arrow key. Keep the down arrow key pressed until the upgrade user software program is launched.



4. The new and current software versions are displayed. Press **Yes** to upgrade to the new software.

Software Upgrade Confirmation

Installation software found on USB mass storage !


Installation software version: X.XX
Current software version: X.XX

Do you want to upgrade it ?

No

Yes

5. An installation progress bar is displayed.

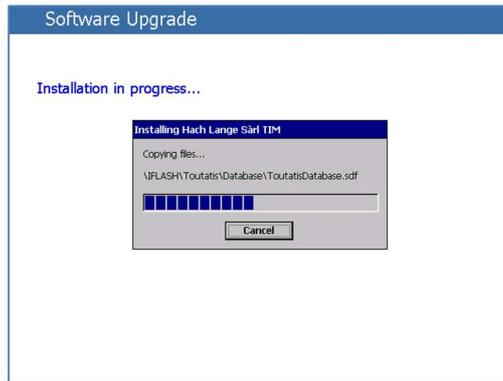

6. After a successful installation, the system will reboot and the new user software will be launched automatically.

